



Examination Light

Outpatient II

The perfect blend of minimal heat output and high color temperature, low-maintenance and top performance.

More physicians choose the Outpatient II than any other light for focusable flood-type light, maximum intensity and shadow-free operation. With minimal heat output and high color temperature, the Outpatient II is the perfect low maintenance, top performing light for nearly every application in the hospital or physician's office.

- Focuses with central removable, autoclavable SteriHandle
- Friction knob allows user to adjust friction on the self-balancing arm positioning to desired tension without a tool
- Cool operation ensured with heat filters and dichroic coatings
- Shadow-free operation provided by three 50-watt halogen bulbs, each with 2000-hour average bulb life
- Bulbs are wired in parallel to ensure continuous operation
- Design allows removal or replacement of SteriHandle with one hand (for ensuring sterile procedures)






burton®



Summary of technical data

OUTPATIENT II

Illuminance:	86,000 Lux at 24" (61 cm)
Color Temperature:	3300 Kelvin
CRI (Color Rendering Index):	94
Light Source:	Halogen
Rated Life:	2000-hour average bulb life
Power:	50 watts
Power:	Floorstand, Single Ceiling, Double Ceiling, Wall, Single Fastrac, Double Fastrac
Product Warranty:	5 Year Limited Warranty

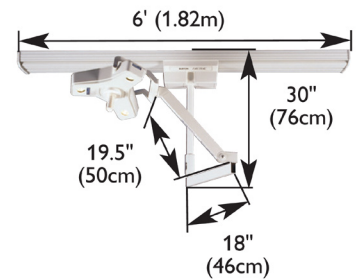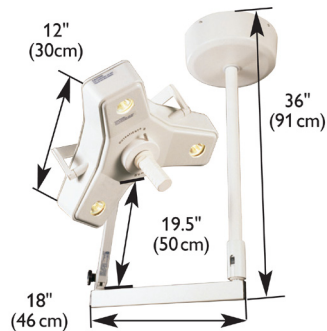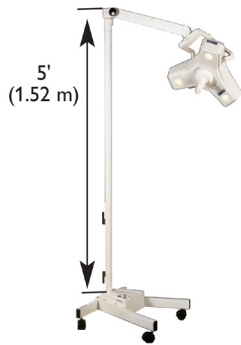

Ordering Information

OUTPATIENT II

120v Model	230/240v Model	230/240v SP Model*	100v Model	Description
OP216FL	OP225FL	O202FL25	OP215FL	Floorstand
OP216SC	OP225SC	O202SC25	OP215SC	Single Ceiling Mount
OP216DC	OP225DC	O202DC25	OP215DC	Double Ceiling Mount
OP216W	OP225W	O302W25	OP215W	Wall Mount
OP216ST	OP225ST		OP215ST	Single Fastrac Mount
OP216DT	OP225DT		OP215DT	Double Fastrac Mount

*SP Model refers to Spanish Language Manuals.

ACCESSORIES

Model No.	Description
0006130PK	Outpatient II Replacement Bulbs, 3/box
1008465	Outpatient II Handle Replacement
0008100PK	Light Handle Covers, Sterile, 25/box
1008259	42" Down Tube
1008060	Special length, max. 58" down tube



© 2022 Burton Medical. All rights reserved. Burton Medical reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

burton®

Burton Medical, LLC
 2300 W Windsor Ct., Suite C
 Addison, IL 60101 USA
 Phone (800) 444-9909, Fax (800) 765-1770
 Burtonmedical.com