



Camera

SmartCam

Making a Practitioners Case Documentation More Efficient.

An autofocus 4K camera for treatment/overhead lights. Burton's SmartCam is specifically designed for lighting mounting. The fast focus and pristine image quality with a Sony IMX415 sensor allows for even easier operation. The SmartCam is a great choice for high level documentation, communication, and live streaming for communication and educational cases.



burton®

Features

- Take Pictures Hands Free and with no Interruption to your Procedure
- Stable Images and Easy to Use (Portable)
- 4K Autofocus Camera


Software

The recording software offers an enhanced and seamless workflow for recording non-obstructed points-of-view videos and still images that can be shown live on a screen in the OR to enhance the documentation of procedures and communication.

With our solution the operatory point of view can be easily shared with the operating team, giving them firsthand sight of the operation. The cameras are USB-connected to ensure continuous operation and recording for the whole procedure and connects to normal medical PCs.

Please visit www.burtonmedical.com/product/novocam/ to download your software.

Burton SmartCam Product Dimensions

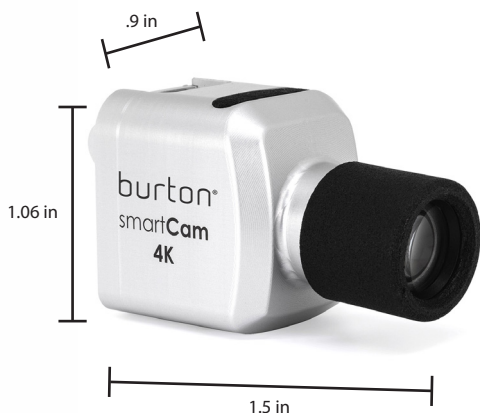

Summary of Technical Data

BURTON SMARTCAM

Sensor:	Sony IMX 415 Starvis
Sensor Resolution:	8.3 MP (4k: 3840x2160)
Sensor Size:	1/2.8"
Video Resolution:	4K, 1080p, 720p; 30 fps
Still Image Resolution:	8.3 MPix
Housing:	Anodized Aluminum
Weight (Body):	1.30
Camera Protocol:	UVC, USB2
Connection:	Micro USB Female on Camera
Lense Options:	25 mm XL NovoCam UHD Lense
Focusing:	Automatic and Manual (SW)
Mounting:	On Light Head
Software Capatibility:	Apple and Windows

Ordering Information

BURTON SMARTCAM

Part Number	Description
SBC4K	Burton SmartCam with Installation Kit Including: USB Foot Control, USB Camera Cable (5m/16ft), USB Extension Cables (2pc 5m/16ft), Cable Clips (30pcs), and Adhesive Ball-Joint for Light Mounting
4000693	Replacement of Cable Clips (30pcs) and Adhesive Ball-Joint for Light Mounting



© 2023 Burton Medical. All rights reserved. Burton Medical reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

burton[®]

Burton Medical, LLC
2300 W Windsor Ct., Suite C
Addison, IL 60101 USA
Phone (800) 444-9909, Fax (800) 765-1770
Burtonmedical.com