



Task Light

Gleamer[®]

Burton's all-purpose halogen light.

If your facility is searching for an all-purpose, rugged, compact halogen light, Burton's Gleamer is your solution. The Gleamer[®] was designed to meet healthcare facilities' most demanding specifications for high-intensity, cool operation, compact design and value. Offering 45,000 lux, a sturdy gooseneck, plus optional 16" (41 cm) extension arm, and superior quality components, the Gleamer provides high-intensity, maximum flexibility and extraordinary reliability. The Gleamer also features a unique venting system, energy-absorbing infrared exchanger and reflector coatings, which filter out infrared for a nearly heat-free beam in a cool, yet compact housing.




burton[®]

Features

- Easy bulb replacement
- Cool housing temperature and cool light beam; light can be in continuous use for hours without causing discomfort to physician or patient
- Convenient positioning with ergonomic positioning ridges
- High-quality, space-saving 4-caster base – safe and flexible portable movement with locking feature in 2 of 4 casters

Gleamer product dimensions

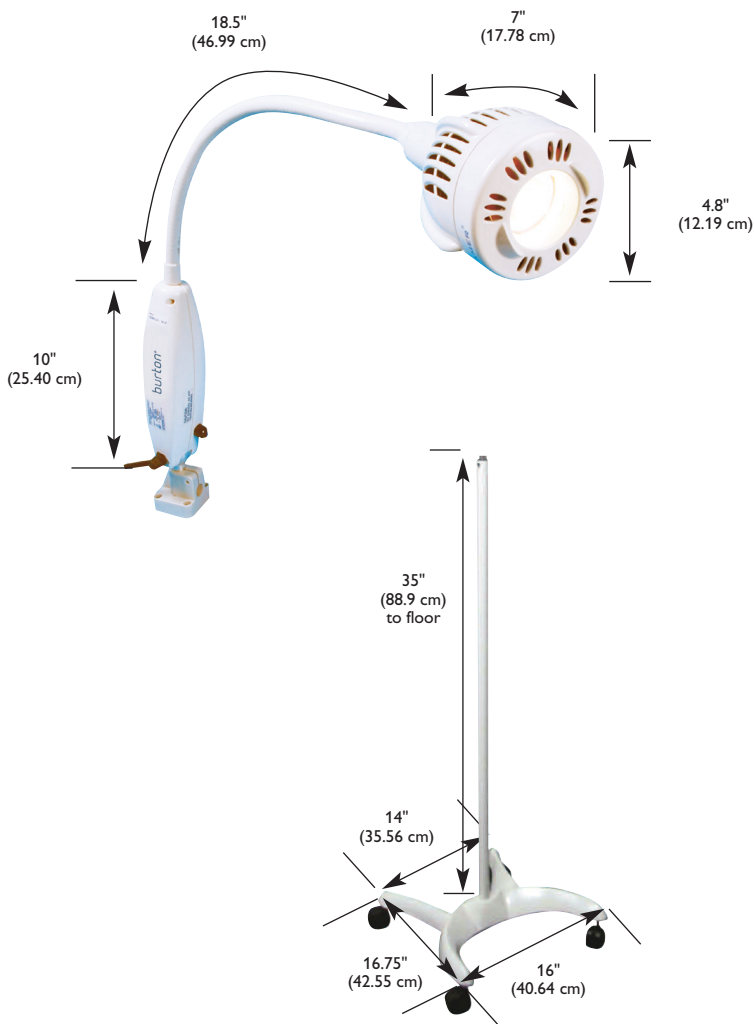

Summary of technical data

GLEAMER

Illuminance:	45,000 Lux (4166 fc) at 18" (46 cm)
Color Temperature:	3700 Kelvin
Light Source:	Halogen
Rated Life:	2,000-hour average bulb life
Power:	35 Watts
Power Supply:	120v model standard; 230/240v models available
Certifications/Approvals:	UL 60601-1/IEC 60601-2-41 Certified
Product Warranty:	5-year limited warranty

Ordering Information

GLEAMER

120v Model	Description
GL30FL	Standard Spot, Floor Stand
GL30W	Standard Spot, Extension Arm & Universal Wall Mount
GLW30FL	Wide Beam Spot, Floor Stand
GLW30W	Wide Beam Spot, Extension Arm & Universal Wall Mount
230/240v Model	Description
GL32FL	Standard Spot, Floor Stand
GL32W	Standard Spot, Extension Arm & Universal Wall Mount
GLW32FL	Wide Beam Spot, Floor Stand
GLW32W	Wide Beam Spot, Extension Arm & Universal Wall Mount

ACCESSORIES

Model	Description
0009600PK	Standard Spot Replacement Bulbs, 4/Box
0009608PK	Wide Beam Spot Replacement Bulbs, 4/Box



© 2018 Burton Medical. All rights reserved. Burton Medical reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

burton®

Burton Medical, LLC
2300 W Windsor Ct., Suite C
Addison, IL 60101 USA
Phone (800) 444-9909, Fax (800) 765-1770
Burtonmedical.com